

M^oDrain[®] Centesis Catheter

The **M^oDrain[®] Centesis Catheter** is an over the needle catheter designed for small abscess fluid aspirations and short-term drainage procedures.


Drainage is accomplished by puncturing the body cavity or organ with a hollow needle and aspirating the contents with a syringe or into an appropriate collection container.

Features:

- **Radiopaque and Kink-Resistant** catheter material
- **Two large oval “skived” holes** to aid in maximum drainage flow and reduction of clogging
- **Unique tapering** provides smooth transition between needle and catheter for easy and precise placement
- **Coated catheter** for easy insertion
- **Echogenic needle tip** to aid in ultrasound visibility
- **Clear needle hub** allows for visual flashback confirmation
- **Porous vent plug** allows air to exit catheter during insertion and protects clinicians from fluid leakage



Smooth needle/catheter transition

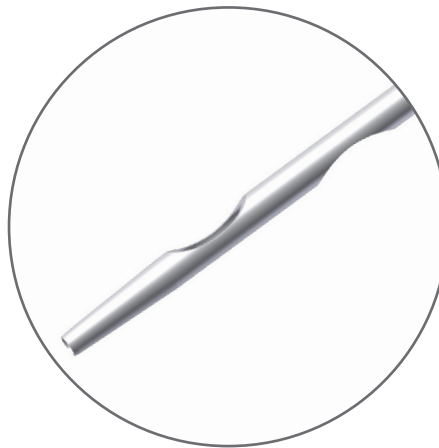

Two large oval “skived” holes



Echogenic tip

M^oDrain[®] Centesis Catheter

Item no.	Size	Qty/box
D40104007	4F x 7 CM Centesis Catheter	10
D40104010	4F x 10 CM Centesis Catheter	10
D40104015	4F x 15 CM Centesis Catheter	10
D40105007	5F x 7 CM Centesis Catheter	10
D40105010	5F x 10 CM Centesis Catheter	10
D40105015	5F x 15 CM Centesis Catheter	10
D40105020	5F x 20 CM Centesis Catheter	10

Contact/Place order: Mermaid Medical A/S · Frydensbergvej 25 · DK-3660 Stenløse · www.mermaidmedical.com · customer@mermaidmedical.com