

Don't trap and drag the clot
when you can **Capture & Retain**



VESALIO is pleased to introduce
FIRST PASS
CLOT REMOVAL TECHNOLOGY

NevaTM

CONTACT

vesalio.com
info@vesalio.com
615.206.7788

PN LV-MKT-002 REV B

ENGINEERED FOR **First Pass Success by Capturing & Retaining All Clot Types**


OPTIMIZED RADIAL FORCE – bi-directional radial force is optimized for clot capture

Proven higher expansive radial force = improved clot integration
 Comparable compressive radial force = safety

DROP ZONE TECHNOLOGY – captures the clot inside the device

Specific areas that allow the clot to drop inside the device, all types of thrombi (soft, organized, hard)

Hybrid stent design with large offset openings optimize clot capture

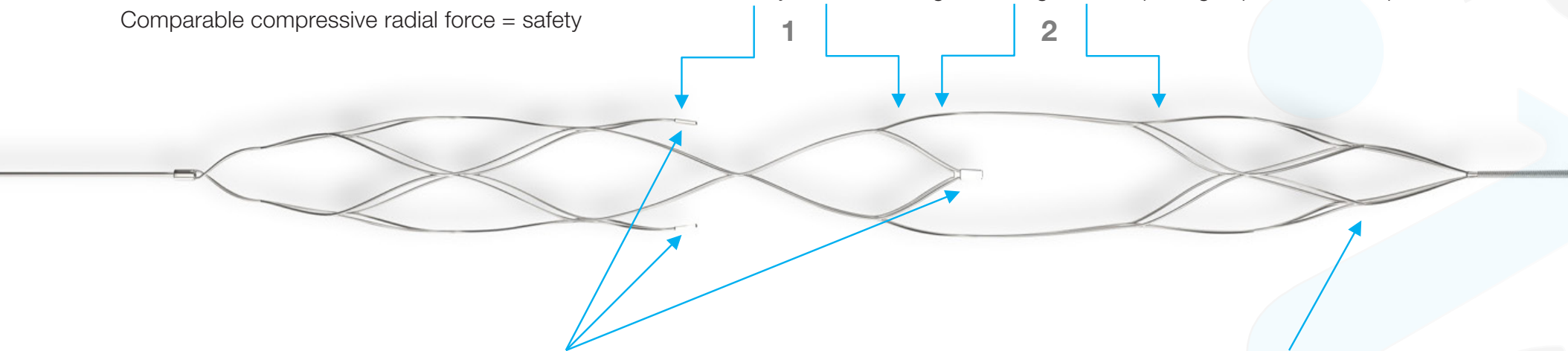

SMART MARKERS – provide physician real time interactive feedback

Placed at the leading edge of the Drop Zones, the tandem markers are offset 90 degrees

Assists operator in considering type of clot, integration and timing & distance of pull

Binary code gives input as to location of clot relative to device

CLOSED DISTAL END – retains the clot and fragments for extraction

Creates an end cap to support the Drop Zones

Reduces downstream migration

PRODUCT FAMILY

Product No.	Product Name	GTIN	Diameter	Device Total Length	Device Working Length	Recommended Vessel Diameter	Pusher Wire Length
30040V-T	NeVa T	00851279008002	4.5 mm	57 mm	37 mm	≤ 4.5 and ≥ 2.0 mm	180.0 cm
30010V-M1	NeVa M1	00851279008019	4.0 mm	48 mm	30 mm	≤ 3.5 and ≥ 2.0 mm	180.0 cm
30020V-MS	NeVa M1-S	00851279008026	4.0 mm	39 mm	22 mm	≤ 3.5 and ≥ 2.0 mm	180.0 cm



Distributed by:

Mermaid Medical Group • Frydensbergvej 25 • DK-3660 Stenløse
 customer@mermaidmedical.com • www.mermaidmedical.com